



"THIS CHAIR HAS MET AND EXCEEDED MY EXPECTATIONS. I'VE GONE MANY MILES AND HOURS ... IN THE ALMOST FOUR YEARS I'VE HAD IT. IT HAS NEVER BROKEN DOWN ON ME OR LET ME DOWN. IT'S ALWAYS GOT ME HOME, THAT'S IMPORTANT. I HAVE CONFIDENCE WHEN I LEAVE ON MY DAILY EXCURSIONS."

Gary Hays - Texas, USA

magic
Mobility™

Frontier V6 All Terrain
**FOR ULTIMATE
OUTDOOR
ADVENTURE**

**FULL LIGHTING
PACKAGES**

Been out having too much fun and it's getting late? Put your lights on. Fitted on rubber brackets.



**UP CLOSE
AND PERSONAL**

Side transfers aren't the only reason you need to flip your armrest out of the way.



REHAB SEATING

Infinite options for infinite body shapes. Our rehab seating range is easily adjustable or can be custom designed to a fixed position. It's whatever you need.

TAKE CONTROL

A full range of advanced control systems and environmental controls are available. Easy to program and easy to use.



 Innovation
In Motion

Exclusive USA distributor since 1999

201 Growth Parkway
Angola, IN 46703
United States

Tel 1-260-665-2769

Email iim@mobility-usa.com

Web mobility-usa.com

V6/USA/DEC19

**BUILT
TO LAST**

**CUSTOMISED FOR
YOUR LIFESTYLE**



**AUSTRALIAN
MADE**

 Innovation
In Motion

Frontier
V6™



FRONTIER V6 ALL TERRAIN

Our most popular powered wheelchair worldwide, the Magic Mobility Frontier V6 All Terrain is ideal for the outdoors, yet its size makes it a capable indoor performer. If performance and versatility is what you are after, the Frontier V6 All Terrain will impress those who prefer the track to the tarmac.

Basic specifications:**

- Length: 41½"
- Width: 28"
- Seat to floor height: 17"
- Turning radius 23"
- Weight limit: 400lbs
340lbs with lift and tilt
- Top speed: 6.2mph
- Motors: 800W Gear in Line
- Wheels: Off-road 14"
(145/70-6) black
- Castors: 2.50-4 black
- Batteries: Group 24, 75 Ah
LA/AGM and Gel (Group 34,
60 Ah Gel option)
- Penny and Giles controller
with LCD display (R-Net)
- Centre post legrest with
one piece footplate
- Sling backrest and rehab
seat pan (from 12"-23" wide)
- WC19 crash test compliant
tie down points.

Rehab positioning:

- Power options: 12" seat lift, tilt
(45 degrees), power recline,
lay down, power elevating legrests,
transfer tilt (8 degrees and up to
2½" floor to seat height reduction),
transfer tilt with tilt only
- Seating and Backrests: MPS,
rehab backrest and adjustable seat
pans and custom seating
- Headrests and mounts: MPS, rehab
- Legrests: Centre post, side mount,
swingaway, power elevating
- Footplates: 1 or 2 piece fixed/
flip up/angle adjustable, any custom
requirement
- A full range of armpads, straps,
pads, harness, restraints, supports
are all available from our extensive
catalogue.

Options and accessories:

- One-click activator
(CJSM 2 or JSM joystick)
- 24V LED Lighting kit
- USB charger
- Jack
- Luggage rack, bag hooks, gun
mounts, bags, pockets, keep safe
secret boxes, tables, drinks holder,
oxy and vent carriers
- Special controls: Attendant, head,
chin, environmental controls all
available
- Hybrid indoor wheel kit: narrower,
grey and black drive wheels and
castors, when off-road capabilities
are not needed. This also reduces
the size of the chair further
- RAM mounting system (tablet and
swing away tablet holder, mobile
phone holder, swing away camera
mount, attendant cup holder, action
camera holder and cup holder).

The right wheelchair for:*

- Adventurous off-road recreation
- Everyday off-road, country living and
working
- Urban outdoor living and working
- Roomy indoor spaces.

Base options:



V6 HYBRID - BLACK WHEELS

"WE'VE BEEN
CREATING CUSTOM
MADE VEHICLES
FOR 25 YEARS, SO
THAT YOU CAN LIVE
THE LIFE YOU
CHOOSE"

Ashley Daff
Director




ATTENTION TO DETAIL

Working together, we can
figure out how to make your
wheelchair uniquely yours.
We know the small touches
make a big difference.

BUILT TO LAST

We never compromise on
quality. With your wheelchair
as your mobility, reliability is
paramount. We'd much rather
spend time improving our
products than repairing them.

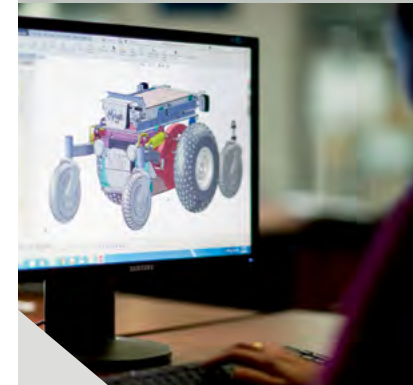

ENGINEERED TO THE LAST BOLT

If we think we can make it better,
we will. Our in house engineering
team are hands on, working
side by side with customers and
therapists. The entire Magic
team work continuously on how
we can improve the experience
from choosing a Magic Mobility
wheelchair to its performance in
the places you want to go.



PROUDLY AUSTRALIAN DESIGNED AND ASSEMBLED IN THE USA

We love that our product is
made in Australia. We love
meeting our customers,
hearing their ever changing
needs and taking it back
to the drawing board.



Magic Mobility +61 3 8791 5600
sales@magicmobility.com.au

Explore with us magicmobility.com.au



Scan for agent locations